

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

For Use in:	Maternity/Women and Children's/ Community Midwifery
By:	Registered Midwives, Senior Maternity Care Assistants, Transitional Care Nursery Nurses
For:	Babies ≥35 weeks gestation and over 24 hours of age being cared for in maternity services
Division responsible for document:	Women and Children
Key words:	Jaundice, Transcutaneous, Neonate
Name of document author:	Debbie Walker, Gayle Richards, Rahul Roy
Job title of document author:	Midwife, Infant feeding co-ordinator, Consultant Neonatologist
Name of document author's Line Manager:	Barbara Jackson, Priya Muthukumar
Job title of author's Line Manager:	Antenatal and Postnatal Matron
Supported by:	Gayle Richards Infant Feeding Coordinator Pam Sizer Community Midwifery Matron Rahul Roy Consultant Neonatologist
Assessed and approved by the:	Clinical Guidelines Assessment Panel (CGAP) If approved by committee or Governance Lead Chair's Action; tick here <input checked="" type="checkbox"/>
Date of approval:	19/06/2020
Ratified by or reported as approved to (if applicable):	Clinical Safety and Effectiveness Sub-Board
To be reviewed before: This document remains current after this date but will be under review	19/06/2023
To be reviewed by:	Debbie Walker, Gayle Richards, Rahul Roy
Reference and / or Trust Docs ID No:	12092
Version No:	3
Compliance links: <i>(is there any NICE related to guidance)</i>	NICE CG 98
If Yes - does the strategy/policy deviate from the recommendations of NICE? If so why?	No

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Version and Document Control:

Version Number	Date of Update	Change Description	Author
2.1	24/04/2020	No clinical changes to document, but due to Covid-19 a short review date has been given to allow for a thorough review of document in the future.	Dr S Mulla and Dr R Roy
2.2	17/06/2020	Thresholds changed. TCB charts have been created using the NICE phototherapy charts by gestational age.	D. Walker RM G. Richards IFC Dr R Roy
3	19/06/2020	Flowcharts tidied , footers amended	Gayle Richards

This is a Controlled Document

Printed copies of this document may not be up to date. Please check the hospital intranet for the latest version and destroy all previous versions.

Contents	Page
Quick Reference Community Midwifery	3
Quick Reference In-patient	4
Objectives	5
Rationale	5
Scope	5
Process	6
Summary of Development	7
References	7
Appendix 1 - Use of JM-105	8
Appendix 2 - Bilirubin (TcB) Threshold Charts	9
Appendix 3 - Daily Operational Check record sheet	11
Appendix 4 - Training Register	14

Definitions of Terms Used / Glossary

Transcutaneous Bilirubinometer	TcB
Serum bilirubin level	SBR
Children's Assessment Unit	CAU
Advanced Neonatal Nurse Practitioner	ANNP
Senior House Officer	SHO
Dräger Jaundice Meter	JM-105

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Community Midwifery

Quick Reference:

Ensure daily check of Dräger Jaundice Meter (JM-105) has been done, passed and recorded before proceeding with test.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Postnatal Inpatient

Quick Reference

Ensure daily check of JM-105 has been done, passed and recorded before proceeding with test.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Objective of Guideline

Aim: To guide the usage of transcutaneous bilirubinometer (TcB) JM-105 as a screening tool for measuring the bilirubin levels in visibly jaundiced well babies ≥ 35 weeks gestation to determine the need for a laboratory serum bilirubin (SBR).

- a. Reduce the need for invasive and painful blood sampling for determining bilirubin level.
- b. Reduce parental anxiety with invasive blood sampling in well babies.
- c. Reduce waiting time for parents.
- d. Aid in referral to CAU as appropriate
- e. Reduce the inconvenience to families arising from the potentially avoidable return visits to hospital.

Rationale for the recommendations

Approximately 60% of term and 80% of preterm babies develop jaundice in the first week of life, and about 10% of breastfed babies are still jaundiced at 1 month of age. In most babies early jaundice is harmless. However, a few babies will develop very high levels of bilirubin, which can be harmful if not treated.

In order to determine the level of jaundice in babies, bilirubin level is currently checked by doctors, ANNP's, midwives and nurses by taking blood via heel prick or venepuncture and then sending it to the laboratory. The transcutaneous bilirubinometer (TcB) is a device, which measures the yellowness of the subcutaneous tissue by a non-invasive method. This method has been recommended by the national guidelines¹ and followed by several regional units in UK. It is a very useful tool to estimate the level of bilirubin in babies at risk or just as part of screening. However, the gold standard method still remains laboratory estimation of serum bilirubin level.²

What is the benefit of using TcB?

- It significantly reduces the number of invasive blood sampling thus reducing workload of clinical staff and laboratory and numbers of babies/families needing to return to hospital.
- More acceptable, reduces anxiety and waiting times for parent.
- Accurate screening tool for visible jaundice.

Scope of practice

- Midwives, ANNPs, Paediatric & Neonatal SHOs and Registrars, Neonatal Nurses, Nursery Nurses and Senior Maternity Care Assistants trained to use JM-105.

Limitations

- The availability of TcBs in the community setting.
- Workload of community staff.
- Availability of staff trained to use the JM-105.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Training

Training will be delivered by designated cascade trainers to those practitioners who will be using the TcB. A training register (Appendix 4) will then be sent to the POCT team who will maintain a register of those staff trained to use the JM-105. It is the practitioner's individual responsibility to ensure they are trained prior to using the JM-105. Training refreshers will be offered on an annual basis by the Practice Development Midwives supported by the Infant Feeding Team.

Process

Babies should be visibly inspected for jaundice at each midwifery contact, as part of the baby examination, ideally naked and in natural light. Inspecting the sclera, gums and blanched skin can aid in the detection of jaundice. Ensure daily check has been done, passed and recorded before using device.

Which babies should be tested with a TcB

- Is visibly jaundiced and more than or equal to 35 weeks gestation
- Is more than 24hrs of age and
- Less than 14 days of age.
- Has not had phototherapy.
- Has jaundice that is worsening, noted either by parent/carer or in professional's opinion.

Which babies should not be tested with a TcB and will need Serum Bilirubin check

- Any baby born under 35 weeks gestation.
- Any baby who required phototherapy
- Any baby less than 24 hours old.
- Any baby over 14 days old.
- Baby who is clinically unwell (e.g. not feeding well, looks septic, dehydrated)

Once a baby has been identified as suitable for TcB reading explain the rationale and procedure to the parents and gain consent to perform the test (See Appendix 1)

After taking the TcB reading (Community)

- Any reading over 250 micromol/litre needs referral to CAU for SBR.
- Take the **highest** of the 3 readings from the JM-105 and compare to the Bilirubin Threshold Chart (Appendix 2).
- If the result is in the green zone continue midwifery-led care taking into account any feeding issues or other red flags.
- If the result is within the yellow zone arrange a repeat in 4 hours. If the result remains yellow refer to CAU for review and SBR, if the result is then green, continue normal care.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

- If the result is in the Red zone refer to CAU for review and SBR – prepare the parents that this could mean a readmission for phototherapy.

After taking the TcB reading (In-Patient)

- Any reading over 250 micromol/litre needs referral for SBR – bleep 0417 (0900-1700) or 0416 (1700-0900) and refer to NICU Assessment Room.
- Take the **highest** reading from the JM-105 and compare to the Bilirubin Threshold Chart (Appendix 2).
- If the result is in the green zone continue midwifery-led care taking into account any feeding issues or other red flags.
- If the result is within the yellow zone repeat TcB in 4 hours. If the result remains yellow refer to on call neonatologist for review and SBR, if the result is then green, continue normal care.
- If the result is in the red zone needs referral to on call neonatologist – bleep 0417 (0900-1700) or 0416 (1700-0900) and refer to NICU Assessment Room.

Summary of development and consultation process undertaken before registration and dissemination

The author listed above drafted this document on behalf of The Infant Feeding Team and Community Midwifery Matron who have agreed the final content. During its development it has been circulated for comment to: Consultant Neonatologists, Midwifery Managers, Infant Feeding Co-ordinators and Postnatal Ward Team Leader and Practice Development Midwives.

1. Bilirubin charts edited to make clear which number were in the red zone > symbols used instead of +.
2. Discrepancies in cleaning clarified.

References

1. NICE clinical guideline 98. Jaundice in Newborn Babies under 28 days. Issue date: October 2016
2. NNUH Clinical Guideline for Non- Invasive Bilirubin measurements using a Transcutaneous Bilirubinometer. [Trustdocs ID No: 12092](#) Issue date February 2017

Associated Documentation

NICE Clinical Guideline 98 Jaundice in Newborn Babies under 28 days


Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Appendix 1


Use of JM-105

How to take a measurement


1. Clean the probe with a Clinell® wipe.

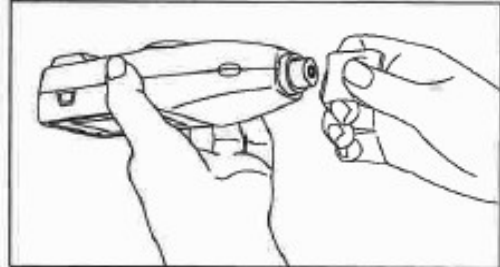

2. Press the power switch on.
3. Select MENU, select MEASURE and press OK.

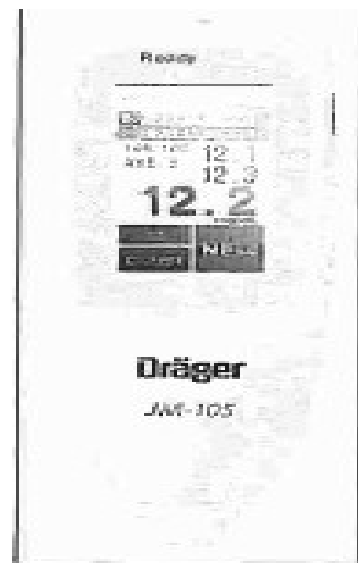

4. Select measurement site: Mid-sternum is preferred but forehead can be used.

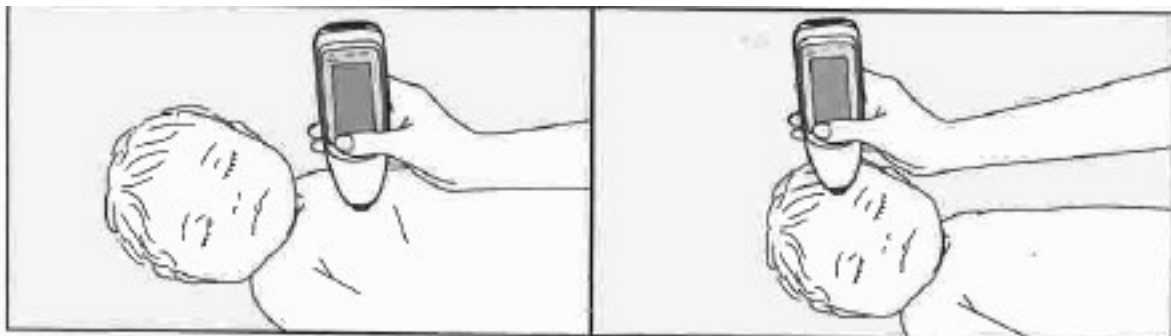

5. Place Jaundice meter JM105 tip flat against the baby's skin, **not at an angle**. Press lightly until you hear a click.
6. Lift the JM-105 from the skin between measurements and pause until green READY light illuminates again. Repeat the testing procedure until the required numbers of measurements (3) have been taken.
7. Record the highest reading on the patient's notes and plot it and record against charts for assessing neonatal jaundice in the neonatal protocol.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

8. Take appropriate action based on the reading (as per the flowchart).

Appendix 2

Bilirubin (TcB) Threshold Charts

These charts have been created using the NICE phototherapy charts by gestational age. Green Zone is greater than 50 below the treatment line, Yellow Zone is within 50 of the treatment line and Red Zone is at or above the treatment line.

35 weeks	Recheck if appears more jaundiced	Requires discussion with oncall CAU paediatrician	Requires urgent referral to CAU for SBR and review
Age in Hours	Bilirubin in micromol/L	Bilirubin in micromol/L	Bilirubin in micromol/L
24 (1 day)	<60	60-110	>110
30	<75	75-125	>125
36	<95	95-145	>145
42	<110	110-160	>160
48 (2 days)	<130	130-180	>180
54	<145	145-195	>195
60	<165	165-215	>215
66	<180	180-230	>230
72 (3 days+)	<200	200-250	>250

36 weeks	Recheck if appears more jaundiced	Requires discussion with oncall CAU paediatrician	Requires urgent referral to CAU for SBR and review
Age in Hours	Bilirubin in micromol/L	Bilirubin in micromol/L	Bilirubin in micromol/L
24 (1 day)	<65	65-115	>115
30	<80	80-130	>130
36	<100	100-150	>150
42	<120	120-170	>170
48 (2 days)	<135	135-185	>185
54	<150	150-200	>200
60	<170	170-220	>220
66	<190	190-240	>240
72 (3 days+)	<210	210-250	>250

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

37 weeks	Recheck if appears more jaundiced	Requires discussion with oncall CAU paediatrician	Requires urgent referral to CAU for SBR and review
Age in Hours	Bilirubin in micromol/L	Bilirubin in micromol/L	Bilirubin in micromol/L
24 (1 day)	<65	65-115	>115
30	<85	85-135	>135
36	<105	105-155	>155
42	<125	125-175	>175
48 (2 days)	<140	140-190	>190
54	<165	165-215	>215
60	<180	180-230	>230
66	<195	195-245	>245
72 (3 days+)	<220	220-250	>250

38+ weeks	Recheck if appears more jaundiced	Requires discussion with oncall CAU paediatrician	Requires urgent referral to CAU for SBR and review
Age in Hours	Bilirubin in micromol/L	Bilirubin in micromol/L	Bilirubin in micromol/L
24 (1 day)	<150	150-200	>200
30	<160	160-210	>210
36	<170	170-220	>220
42	<190	190-240	>240
48 (2 days)	<200	200-250	>250
54	<210	210-250	>250
60	<225	225-250	>250
66	<235	235-250	>250
72 (3 days+)	<250 (Alert and feeding effectively)	<250 (Lethargic or Feeding Issues)	>250

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Appendix 3

Drager JM-105 Jaundice Meter: Daily Operational Check Record Sheet

Month:

SN:

Date	L value range: -2.9 to -0.9	S value range: -2.6 to -0.6	Δ value range: -0.8 to +0.2	?within range	Signature
1					
2					
3					
4					
5					
6					
7					
8					
9					
10					
11					
12					
13					
14					
15					
16					
17					
18					
19					
20					
21					
22					
23					
24					
25					
26					
27					
28					
29					
30					
31					

L value = measured value of long optical path

S value = measured value of short optical path

Δ value = difference between L and S values

- If any value is out of range, clean the checker and probe. Repeat the measurement.
- If values are still or repeatedly out of range, do not use the device and contact Point Of Care Testing Department (POCT) on extension 2969 or 01603 286969 extension 2710 or 01603 287710

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

Preparation of the JM-105 for Use:

Perform the Daily Operational Checkout Procedure

1. Remove the JM-105 from the docking station.
2. Press the power switch on.
3. Select 'CHECKER' and touch 'OK' to save selection.
4. Open the checker lid on the charging unit.
5. When the green READY light illuminates, place the tip of the meter perpendicular on the reading checker circle. Press down until you hear a click.
6. The display screen shows the "L" (long), "S" (short), and Delta values.
7. The meter must read within the reference values posted under the checker lid.
 - a. If so, the meter is ready to use.
 - b. If not, clean the tip and repeat.

If values are still out of range, do **NOT** use the meter and contact the POCT Co-ordinator. extension 2969 or 01603 286969 extension 2710 or 01603 287710

Storing the JM-105

- Clean the measuring probe.
- Place the meter in the docking station.
- Leave the device in the docking station when not in use.
- While the meter is in the docking station, power remains on. After 1 minute the screen goes blank; the power remains on.

Switching off the JM-105

To switch off power when the device is not in the docking station, press the power button and hold it down for 1 second.

Precautions

- **Cleaning**
 - The probe must be cleaned with Clinell wipes prior to testing each baby.
 - In the event of an isolated infectious baby, chlorine (Haztabs 1,000ppm) to be used on device.
 - If device is contaminated with blood or body fluids, chlorine (Haztabs 10,000 ppm) to be used.

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

- A device check must be performed daily. If a check has not been carried out during the current day, the **MEASURE READING CHECKER** message appears for 3 secs when the device is switched on.
 - To clear the message, perform the operational checkout procedure.
- If check values are repeatedly out of range, do not use the JM-105; contact the POCT Co-ordinator on extension 2969 or 01603 286969 extension 2710 or 01603 287710

Limitations

- Only staff who have been trained by the Dräger Training Specialist or Link Trainers are authorised to perform patient testing.
- Use only on infants up to 14 days of age.
- Use only on infants of ≥ 35 weeks gestation at birth.
- Use only the sternum location when taking measurements (or forehead, at hospital sites only), where a sufficient amount of blood is circulated.
 - Clinical studies show consistently better results with measurements taken at the sternum versus the forehead.

Problem solving – error messages

Error code	Cause	Action
ERROR0 1	Measured value is abnormal	Repeat the measurement
ERROR0 3	RAM error	Switch off and do not use
ERROR0 3 ERROR0 4 ERROR0 5	Averaging failure Hardware failure	Switch OFF power Wait 10s Switch ON
ERROR0 4	Memory error	Switch off device and do not use
ERROR0 5	Insufficient charge, circuit error	Charge the device If failure continues, switch off and do not use
ERROR0 6	Calibration data error	Switch off the device and do not use

Clinical Guideline/Standard Operating Procedure for the for the Use of Transcutaneous Bilirubinometers in Neonatal and Maternity Services

